


# BD Stent

SX-ELLA Stent Esophageal Degradable BD



# BD Stent

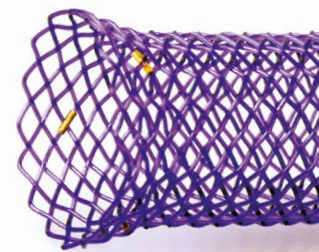

Biodegradation

Unique delivery system

Excellent flexibility

Radiopaque markers



## Basic description

# SX-ELLA Stent Esophageal Degradable BD – BD Stent

BD Stent is designed for relieving benign esophageal lesions, namely:

- stenosis (peptic, anastomotic or caustic) refractory to standard therapy,
- achalasia refractory to standard therapy.

## Features / Benefits





- ✓ Biodegradable material - Polydioxanon / Disintegration of stent usually occurs after 11 – 12 weeks following implantation / No need for removal
- ✓ PULL-back delivery system / Splittable olive tip / Automatic detachment while stent starts releasing / Prevention of olive entrapment / Easy insertion and safe removal of delivery system
- ✓ Radiopaque markers / Accurate stent positioning
- ✓ Anti-migration design / Dual flared design / Reduced migration rates
- ✓ Integrity and radial force of stent is maintained for 6 – 8 weeks following implantation / Extended period of dilation compared to conventional methods
- ✓ Simple compression procedure / Compression tool available

## Additional description

Stent is supposed to be compressed by compression tool into delivery system just prior to implantation

6 mm / 9.4 mm (18 F / 28 F) delivery system, available in standard length of 75 cm

Recommended guide wire: 0.035" (0.89 mm) ultra stiff / Minimum length 220 cm

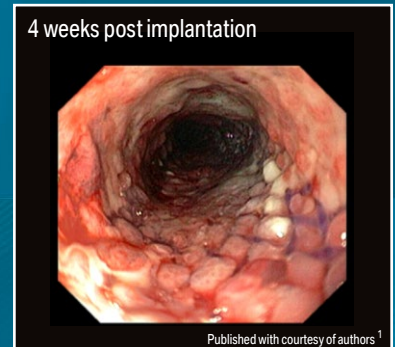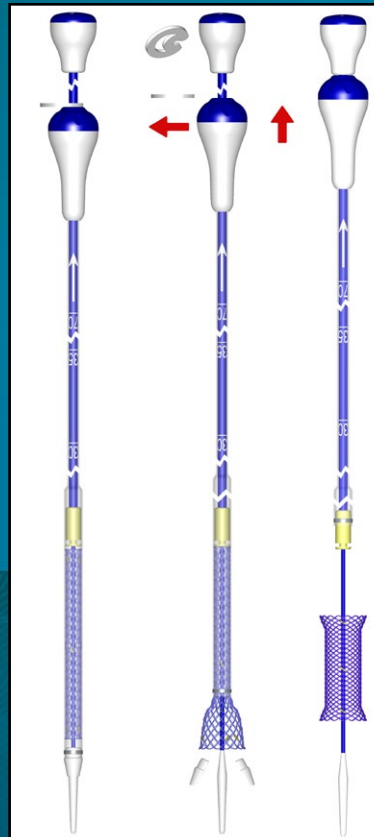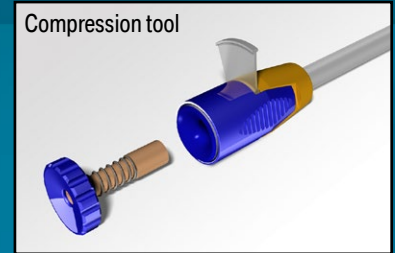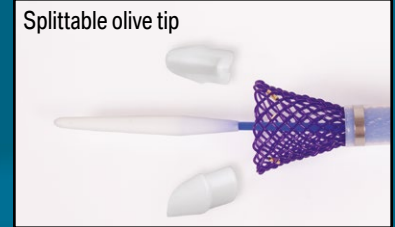

Published with courtesy of authors<sup>1</sup>



## Available sizes

REF No.	Stent body diameter [mm]	Stent flares diameter [mm]	Stent length [mm]
019-10A-23/18/23-060	18	23	60
019-10A-23/18/23-080			80
019-10A-23/18/23-100			100
019-10A-25/20/25-060	20	25	60
019-10A-25/20/25-080			80
019-10A-25/20/25-100			100
019-10A-28/23/28-060	23	28	60
019-10A-28/23/28-080			80
019-10A-28/23/28-100			100
019-10A-31/25/31-060	25	31	60
019-10A-31/25/31-080			80
019-10A-31/25/31-100			100
019-10A-31/25/31-135			135

Distributor:

Manufacturer:

**ELLA-CS, s.r.o.**

Milady Horakove 504/45, Trebes,  
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655

E-mail: info@ellacs.eu URL: www.ellacs.eu