

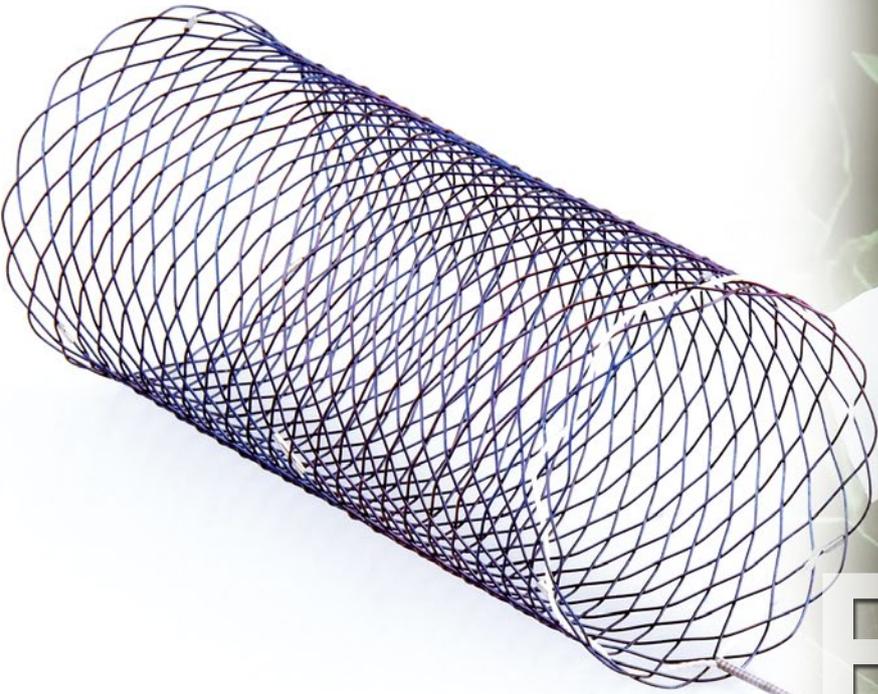


# Enterella

## SX-ELLA Stent Colorectal

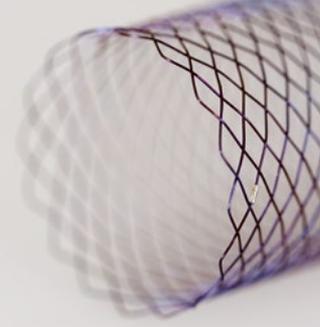

# Enterella

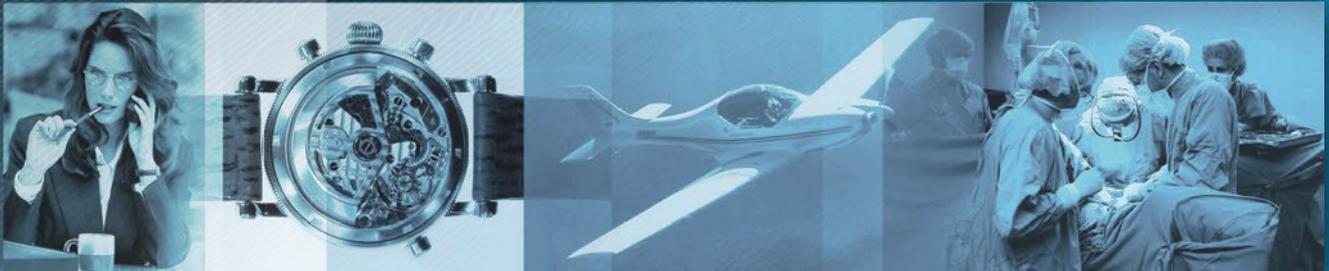

Easy deployment

Excellent flexibility

Atraumatic stent ends

Smooth alignment around flexures



## Basic description

### SX-ELLA Stent Colorectal – Enterella

Stent system for insertion under fluoroscopic or endoscopic guidance. Self-expandable nitinol stent preloaded in a pull-back delivery system.

## Indications

Enterella Colorectal is indicated for palliative treatment of the obstructions or strictures caused by malignant tumors which affect the rectum, the colon sigmoideum or distal part of the descending colon. It may be also used for temporary relieving such obstructions prior to the elective surgery.

## Contraindications

- ☑ Bowel perforation
- ☑ Peritoneal abscess
- ☑ Severe ascites
- ☑ Severe coagulopathy
- ☑ Enteral ischemia
- ☑ Placement of intestinal obstruction / stricture outside the anatomic structures defined in "Indications"

## Features and benefits

- ☑ Atraumatic stent ends - effective prevention of trauma to the bowel wall.
- ☑ Excellent flexibility - precise alignment around anatomical curves.
- ☑ Radiopaque markers at both stent ends and at the midpoint - improved visibility; accurate stent positioning.
- ☑ Plastic repositioning loop at aboral stent end - possibility of stent repositioning within five days after implantation.


## Additional description

15 F (5.08 mm) delivery system available in a standard length of 95 cm (135 cm delivery system is available on request).

Radiopaque, atraumatic olive tip of the delivery system.

Recommended guide wire: 0.035" (0.89 mm) ultra stiff.


Excellent flexibility

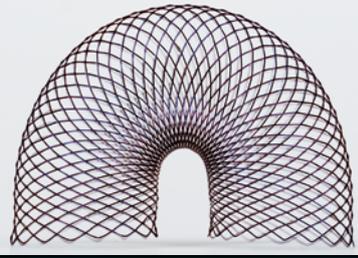


Plastic repositioning loop



Atraumatic stent ends

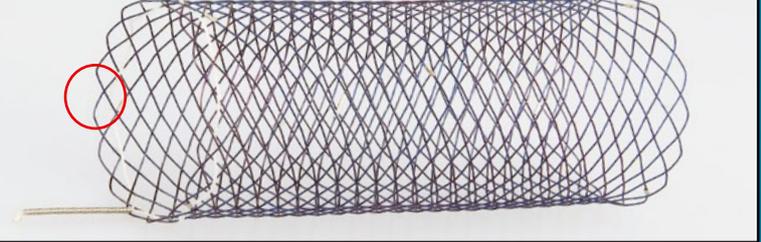


Easy deployment

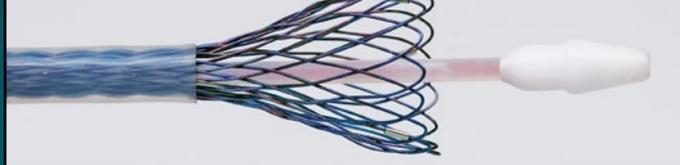

## Available sizes

REF No. Uncovered version	Nominal (unconstrained) diameter [mm]	Nominal (unconstrained) length [mm]
019-05-22-082	22	82
019-05-22-090	22	90
019-05-22-113	22	113
019-05-22-135	22	135
019-05-25-082	25	82
019-05-25-090	25	90
019-05-25-113	25	113
019-05-25-135	25	135
019-05-30-075	30	75
019-05-30-088	30	88
019-05-30-112	30	112
019-05-30-123	30	123
019-05-30-136	30	136

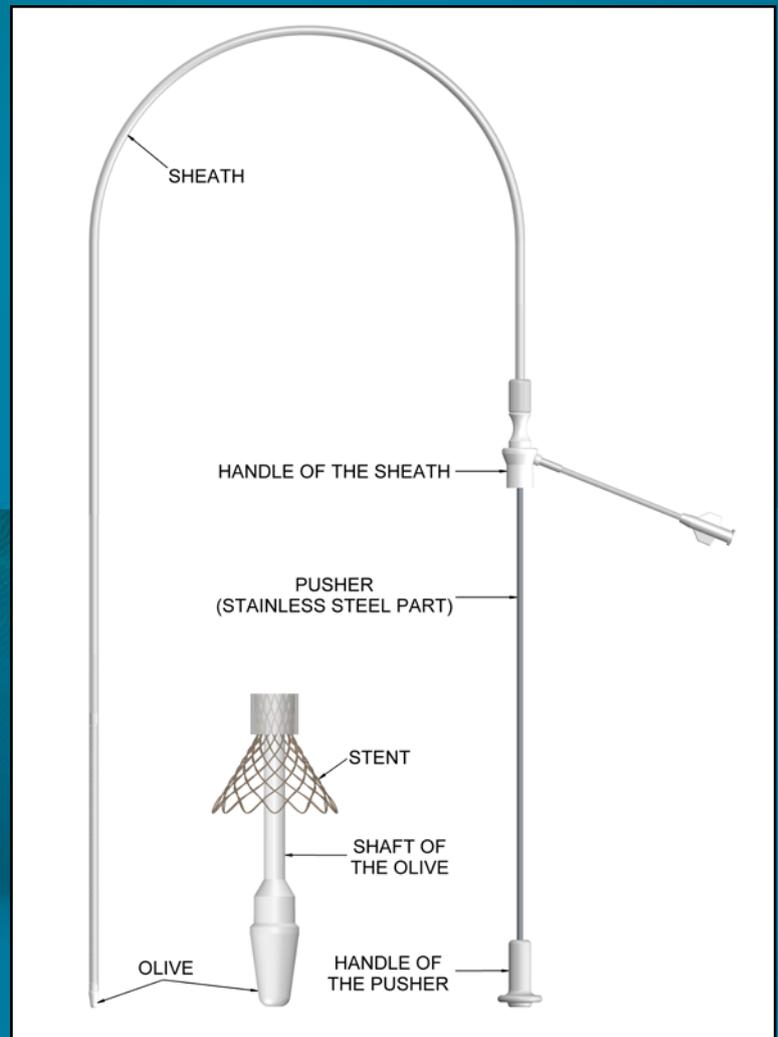

Distributor:

Manufacturer:

**ELLA-CS, s.r.o.**

Milady Horakove 504/45, Trebes  
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655  
E-mail: info@ellacs.eu URL: www.ellacs.eu