



Enterella

SX-ELLA Stent Pyloroduodenal

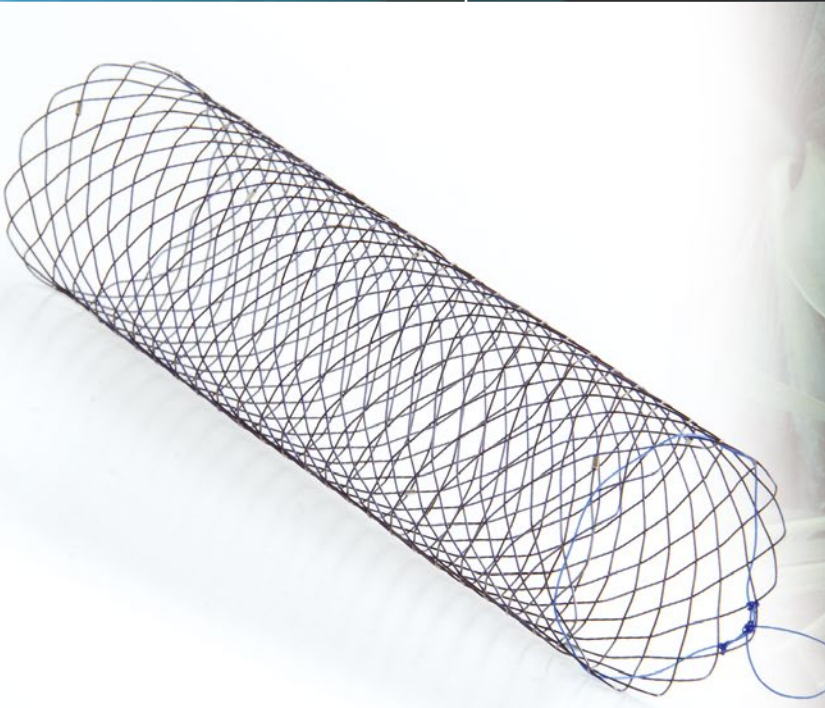

Enterella

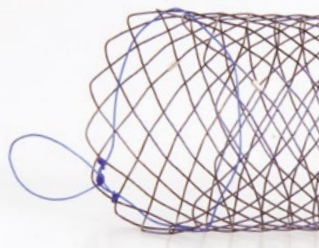

Easy deployment

Excellent flexibility

Atraumatic stent ends

Smooth alignment around flexures



Basic description

SX-ELLA Stent Pyloroduodenal – Enterella

Stent system for fluoroscopic and/or endoscopic insertion.
Self-expandable nitinol stent preloaded in a pull-back delivery system.

Indications

Enterella Pyloroduodenal is indicated for palliative treatment of gastric outlet syndrome or duodenal strictures and/or obstructions caused by malignant tumors.

Contraindications

- ☑ Enteral ischemia
- ☑ Bowel perforation
- ☑ Peritoneal abscess
- ☑ Severe ascites
- ☑ Severe coagulopathy


Features / Benefits

- ☑ Atraumatic stent ends - Effective prevention of trauma to the gastric/duodenal wall.
- ☑ Excellent flexibility - Precise conformation to the anatomical curves in a patient.
- ☑ Radiopaque markers at both stent ends and at the midpoint - Improved visibility; accurate stent positioning.
- ☑ Plastic repositioning loop at the proximal stent end - Possibility of stent repositioning within five days following implantation.
- ☑ Braided kink-resistant delivery system - Superior system introduction and stent deployment.


Additional description

10.5 F (3.5 mm) delivery system available in a standard length of 210 cm allowing passage through the working channel of an endoscope.
Radiopaque, atraumatic olive tip of the delivery system.
Recommended guide wire: 0.035" (0.89 mm).


Kink-resistant delivery system

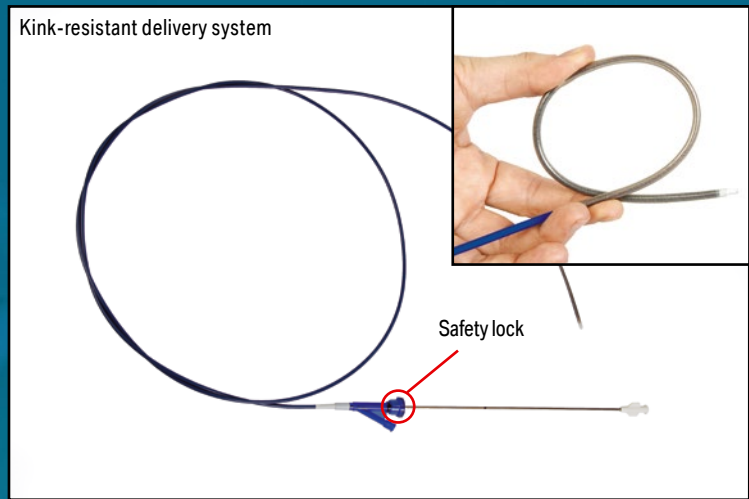


Yellow visual marker and radiopaque ring

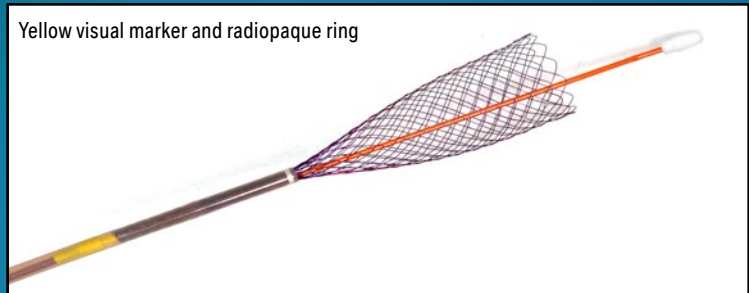


Radiopaque markers

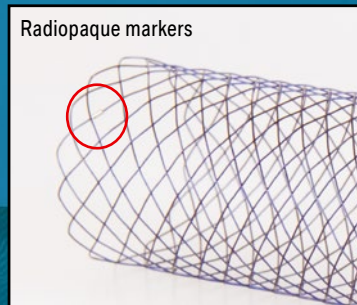


Plastic repositioning loop

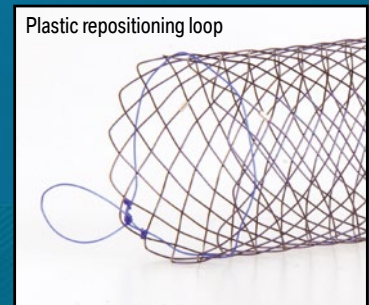

Atraumatic stent ends

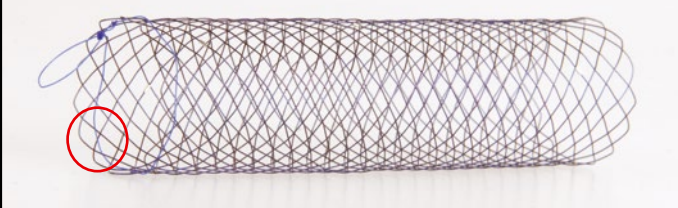

Easy deployment

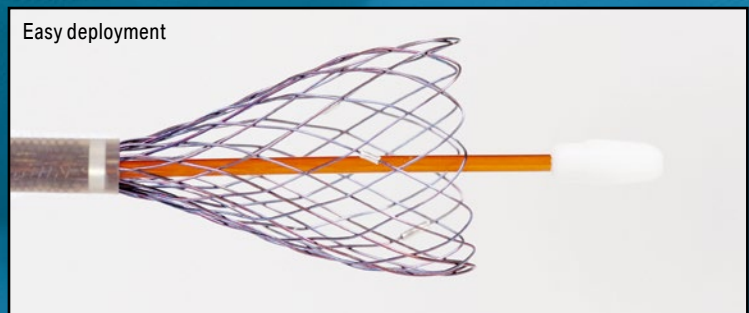

Available sizes

REF No.	Nominal (unconstrained) diameter [mm]	Nominal (unconstrained) length [mm]
019-04-20-082	20	82
019-04-20-090	20	90
019-04-20-113	20	113
019-04-20-135	20	135
019-04-22-082	22	82
019-04-22-090	22	90
019-04-22-113	22	113
019-04-22-135	22	135
019-04-25-082	25	82
019-04-25-090	25	90
019-04-25-113	25	113
019-04-25-135	25	135

Distributor:

Manufacturer:

ELLA-CS, s.r.o.

Milady Horákové 504/45, Třebeš
500 06 Hradec Králové, Czech Republic

Tel.: +420 495 279 111 E-mail: info@ellacs.eu
www.ellacs.eu